



HEAT SEALED POWDER COAT STEEL 16 SCOPE TAMBOUR DOOR DRYING CABINET, KEY LOCK

Scopes not
included

SC8036TKDP-16

SCOPE CABINET FEATURES

- Positive pressure fully redundant fan system
- Includes two high volume intake fans and one exhaust fan
- Fan run time is fully adjustable and activated each time cabinet is closed
- HEPA filter removes 99.7% of particulates, .03 Microns
- Holds 16 scopes with a secure adjustable mounting system (other configurations available)
- Convenient replacement of HEPA filters with easily removed thumb screws
- Medical grade electrical system
- Roll-up tambour door with key lock
- Easy to clean removable drip tray
- Cabinets are delivered fully assembled
- Harloff has Industry Best 12 year limited manufacturers warranty for Sales to US and Canada
- Cabinet dimensions: 93.625"H x 36"W x 24"D

COLORS

Available in more than twenty color options.

View online: <https://www.harloff.com/finish-options/>

The Sure-Dry Scope Cabinet includes the industry best Positive Pressure HEPA Filtered Drying system. Inside air is replaced inside the cabinet four times per minute.



Harloff scope cabinets meet or exceed Joint Commission, AORN, SGNA and AAMI guidelines.



Engineered &
Manufactured
in the U.S.



View similar products online at
<https://www.harloff.com/suredry-scope-cabinets/>

**Toll Free: (800) 433-4064 • medicalsales@harloff.com
www.harloff.com • Ph: (719) 637-0300**

10 Reasons to Choose **HARLOFF** for Your Scope Cabinets



Industry Best Positive Pressure HEPA Filtered Drying System. Designed to add Dri-Scope Aid® now or later with easy “Plug & Play” installation.

1

A Smart Choice to Meet & Exceed Current SGNA, AORN and AAMI Guidelines

Proudly Made in the USA



Harloff products are proudly **Designed, Engineered and Manufactured in the United States** at our Colorado Springs, CO facility.

6


Dri-Scope Aid® Fully Retrofittable in Field

All SureDry scope cabinets are fully designed to be field installed with the industry best Dri-Scope Aid® lumen drying system, which provides HEPA filtered air to each internal scope channel.

2

Easy Access to HEPA Filters

Change HEPA filters fast **using simple thumb screws**. Save time and minimize cycle interruptions.

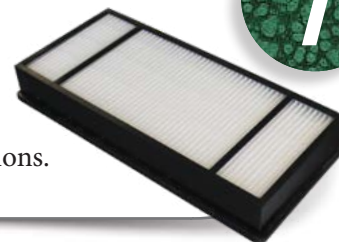

7



LED Lighting Option

Brighten your cabinet and easily inspect scopes with fully dimmable optional overhead LED lighting which triggers on door open/close.

3



Easy Lifting Tambour Door

There's no need to bend down to lift the scope cabinet door. **With an ergonomic handle built right into the middle of the tambour door**, opening and closing your scope cabinet requires no stooping and minimal effort.

8

Hooks and Holders for Specialized Scopes

Major brands to include **Fuji, Olympus, Pentax, Karl Storz and others.**



4

Scope Cabinet Finish Options

Choose Stainless Steel, Heat Sealed Powder Coat in High Gloss or Textured Low Gloss, or upgrade to an **antimicrobial finish**.

9

12 Year Warranty

Our Manufacturer's Warranty is Unmatched in the industry. We build our cabinets to last and we stand behind them.

5



Superior Color Choices

Over 20 Colors to choose from. Cabinets have a heat sealed powder coat finish inside and out. View color options online at www.harloff.com/finish-options/

10